

Application Usage and Features

Launching The Application

In order to Launch the decimal Ordering Tool app you must use the decimal Launcher to authenticate and open your desired application.

NOTE: You will need your .decmlal Direct user Log in.

1. Download and install the decimal Launcher.
2. Log in with your .decmlal credentials.
3. Find decimal Ordering Tool and download/launch the application.

For further instructions please refer to the [decimal Launcher User Guide](#).


Ordering a Device

The decimal Ordering Tool supports selecting and ordering the following file types, part types, and materials:

File Type	Part Type	Material
STL	Bolus	Blue Wax
		FlexiBol (silicone)

Order Information

Once a supported file type is opened the order information can be filled in. See below for a description of each order field and how it will be used for the order:

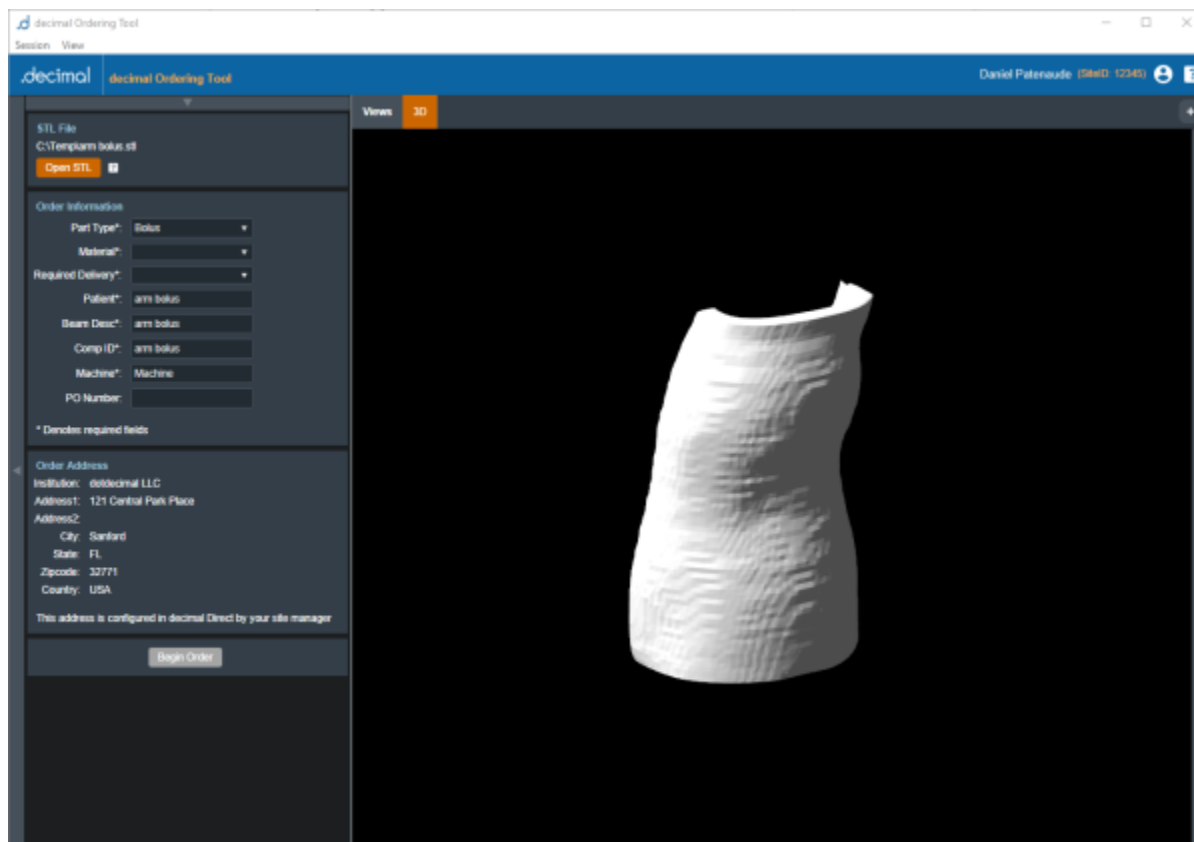

Field	Description
Part Type	Required Field. The specified part type to be ordered
Material	Required Field. The material for the order based on the selected part type
Required Delivery	Required Field. The manufacturing and shipping speed desired: <ul style="list-style-type: none"> • Standard: Overnight delivery by end of day • Rush: Overnight delivery before noon Note: extra charges may apply for Rush delivery Refer to the decimal Direct Shipping Estimate Calculator for a detailed estimate for when your part will ship from .decimal.
Patient	Required Field. An anonymous identifying label for the order's patient. This defaults to the selected file name. Note: .decimal requires this to not contain plain text PHI/PII.
Beam Description	Required Field. A user desired description for this beam/device present on order reports. This defaults to the selected file name. Note: .decimal requires this to not contain plain text PHI/PII.
Comp ID	Required Field. A user desired label for this device that will be included on the hardware labeling and present on order reports. This defaults to the selected file name. Note: .decimal requires this to not contain plain text PHI/PII.


Field	Description
Machine	Required Field. Any text label that describes the machine being used for the part. This defaults to "Machine" and is not used for order reports or order processing.
PO Number	Optional Field. If applicable, an existing Purchase Order number that will be included and referenced on the final Invoice for the order

Order Address

The order address is automatically configured for your site within [decimal Direct](#) by your Site Manager. Please contact your Site Manager or .decimal Customer Service (Phone: 1.800.255.1613, E-Mail: customersupport@dotdecimal.com) if you require changing this address.

Placing an Order

Once all required order fields have been filled in the Begin Order button will enable allowing you to order the device to [decimal Direct](#).

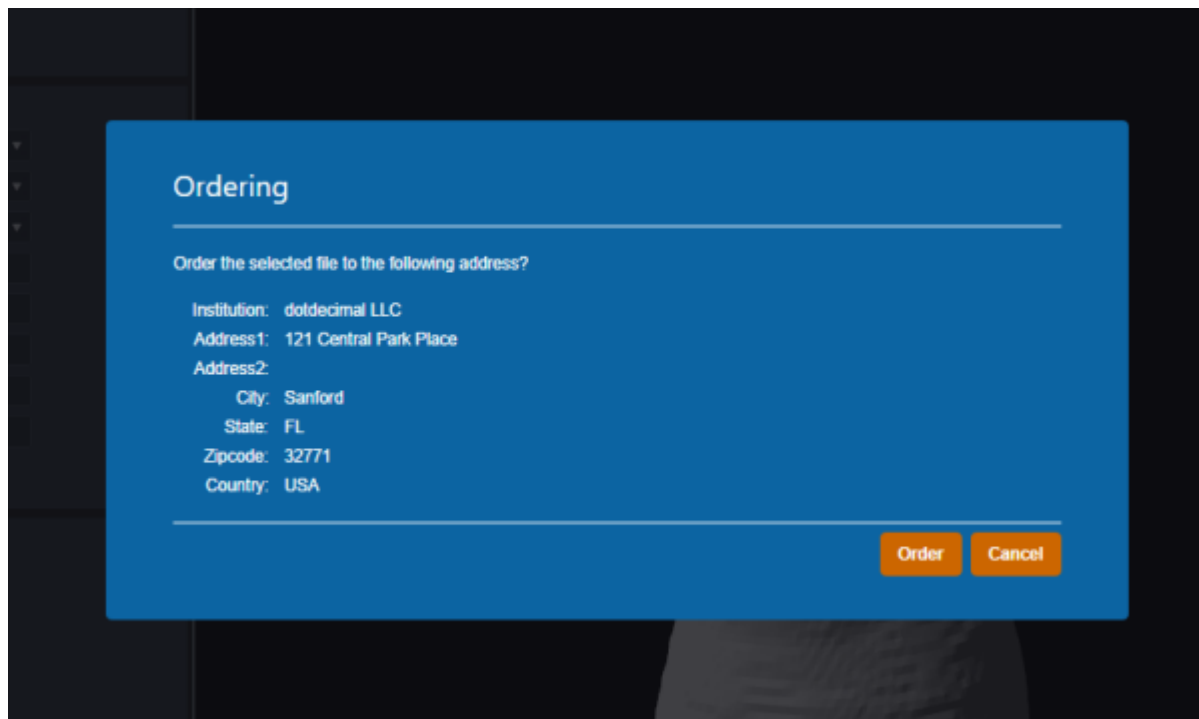A screenshot of a web application showing a modal window titled "Ordering". The modal has a blue header and a white body. Inside the body, it says "Order the selected file to the following address?". Below this, it lists the following information: Institution: dotdecimal LLC, Address1: 121 Central Park Place, Address2: (blank), City: Sanford, State: FL, Zipcode: 32771, and Country: USA. At the bottom right of the modal, there are two orange buttons: "Order" and "Cancel".

Ordering

Order the selected file to the following address?

Institution: dotdecimal LLC
Address1: 121 Central Park Place
Address2:
City: Sanford
State: FL
Zipcode: 32771
Country: USA

[Order](#) [Cancel](#)

Confirm the order address and

Application Settings

The following settings are present within the decimal Ordering Tool settings UI accessed within the View

→ Settings menu:

General

- **Landing Page**
 - **Do not show the app welcome page:** Allows the user to specify whether or not to view the application landing page when decimal eRT is opened.
- **UI Theme**
 - Set the visual theme of the application (Note: the 'blue' theme is the currently supported theme; a 'dark' theme is available, but is a deprecated theme and not routinely updated at this time).

Application

There are no specific application settings for the decimal Ordering Tool at this time.

Cache

The local disk cache is where calculation results are stored on the client computer.

Users can choose to manually clear the disk cache or let it clear oldest data as the cache fills. Refer to [Data Management Storage and Caching](#) for more information.

Security

- **Inactivity Timeout**
 - Specifies the time out period (in minutes) to automatically log the user out of the application. Timeout is based on mouse/keyboard interaction within the user interface. Note: This inactivity timeout is different than the user session token received from the [decimal Launcher](#).
- **SSL Settings**
 - **Disable SSL Certificate Revocation Check:** Enabling this option will disable SSL Certificate Revocation checking for proxy networks. This is not recommended unless your facility network is using a proxy that is causing SSL certificate validation errors within the application.

From:
<http://apps.dotdecimal.com/> - **decimal App Documentation**

Permanent link:
http://apps.dotdecimal.com/doku.php?id=decimal_ordering_tool:userguide:tutorials&rev=1654279228

Last update: **2022/06/03 18:00**

