

# Application Usage and Features

## Launching The Application

In order to Launch the decimal Ordering Tool app you must use the decimal Launcher to authenticate and open your desired application.

**NOTE:** You will need your .decmlal Direct user Log in.

1. Download and install the decimal Launcher.
2. Log in with your .decmlal credentials.
3. Find decimal Ordering Tool and download/launch the application.

For further instructions please refer to the [decimal Launcher User Guide](#).

---

## Ordering a Device


The decimal Ordering Tool supports selecting and ordering the following file types, part types, and materials:

File Type	Part Type	Material
STL*	Bolus	Blue Wax
		FlexiBol (silicone)

\*Note: STL files are required to be binary format and not ascii.

## Order Information

Once a supported file type is opened the order information can be filled in. See below for a description of each order field and how it will be used for the order:

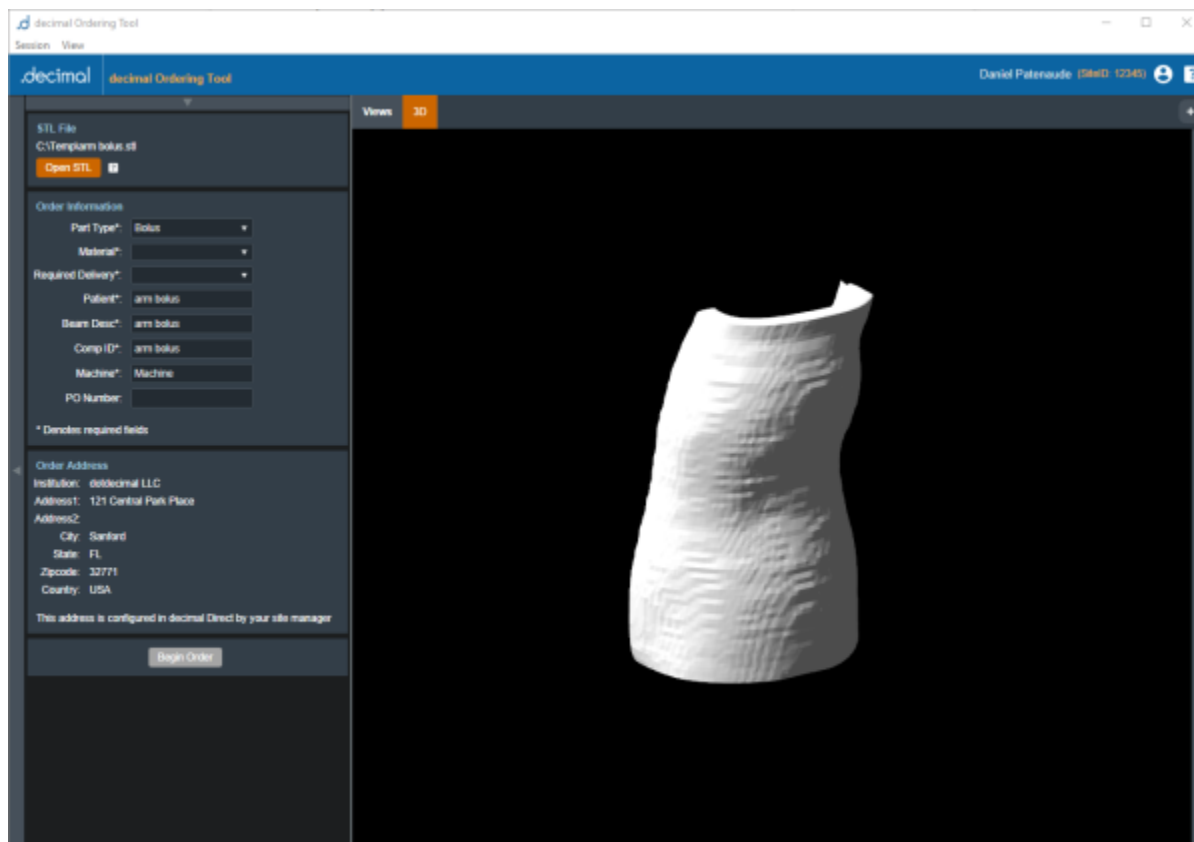

Field	Description
Part Type	Required Field. The specified part type to be ordered
Material	Required Field. The material for the order based on the selected part type
Required Delivery	Required Field. The shipping speed desired*: <ul style="list-style-type: none"> <li>• <b>Standard:</b> Overnight delivery by end of day</li> <li>• <b>Rush:</b> Overnight delivery before noon</li> </ul> Note: extra charges may apply for Rush delivery *Refer to the <a href="#">decimal Direct Shipping Estimate Calculator</a> for a detailed estimate for when your part will ship from .decimal.
Patient	Required Field (limit 50 characters). An anonymous identifying label for the order's patient. This defaults to the selected file name. <b>Note: .decimal requires this to not contain plain text PHI/PII.</b>
Beam Description	Required Field (limit 50 characters). A user desired description for this beam/device present on order reports. This defaults to the selected file name. <b>Note: .decimal requires this to not contain plain text PHI/PII.</b>
Comp ID	Required Field (limit 50 characters). A user desired label for this device that will be included on the hardware labeling and present on order reports. This defaults to the selected file name. <b>Note: .decimal requires this to not contain plain text PHI/PII.</b>

Field	Description
Machine	Required Field. Any text label that describes the machine being used for the part. This defaults to "Machine" and is not used for order reports or order processing.
PO Number	Optional Field. If applicable, an existing Purchase Order number that will be included and referenced on the final Invoice for the order

**PHI/PII Warning**

Order information fields can ship up on invoices, packing slips, and other places not guaranteed to be protected by authorized access. As such, these fields should not contain any PHI/PII.

## Order Address

The order address is automatically configured for your site within [decimal Direct](#) by your Site Manager. Please contact your Site Manager or .decimal Customer Service (Phone: 1.800.255.1613, E-Mail: [customersupport@dotdecimal.com](mailto:customersupport@dotdecimal.com)) if you require changing this address.

## Placing an Order

Once all required order fields have been filled in the Begin Order button will enable allowing you to order the device to [decimal Direct](#).

**Ordering**


---

Order the selected file to the following address?

Institution: dotdecimal LLC  
Address1: 121 Central Park Place  
Address2:  
City: Sanford  
State: FL  
Zipcode: 32771  
Country: USA


[Order](#) [Cancel](#)

Confirm the order address and place the order and the order will be uploaded, validated, and placed to [decimal Direct](#). Once the order is successfully placed you will receive a successful order confirmation message within the application.

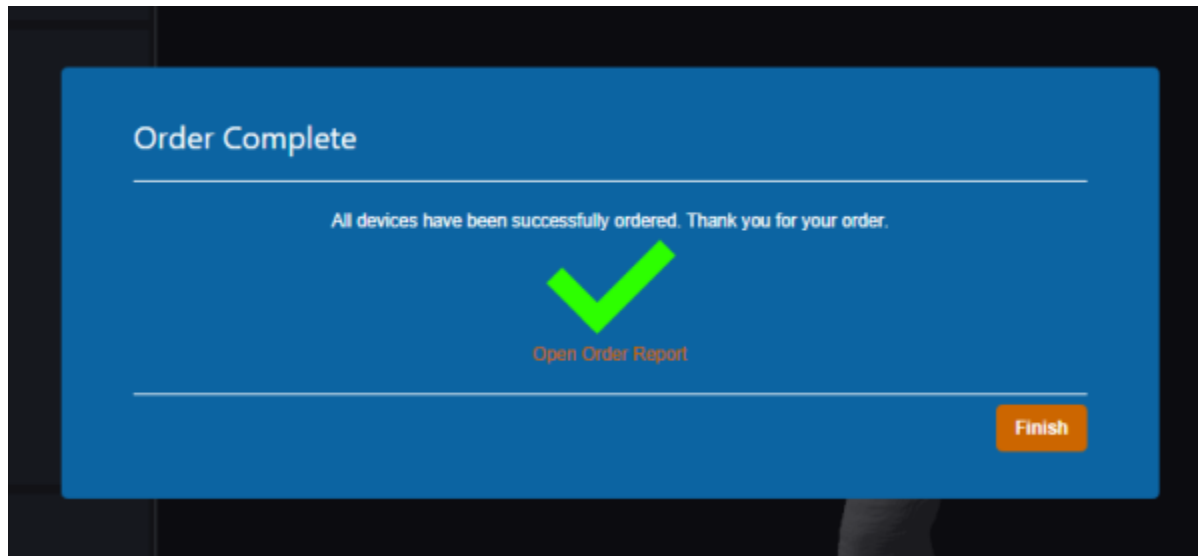

## Ordering Errors

If there was an error placing your order the application will display the error, details of the error, and provide a link to the [decimal Direct ordering page](#). Clicking this link will automatically navigate you and log you in to decimal Direct.

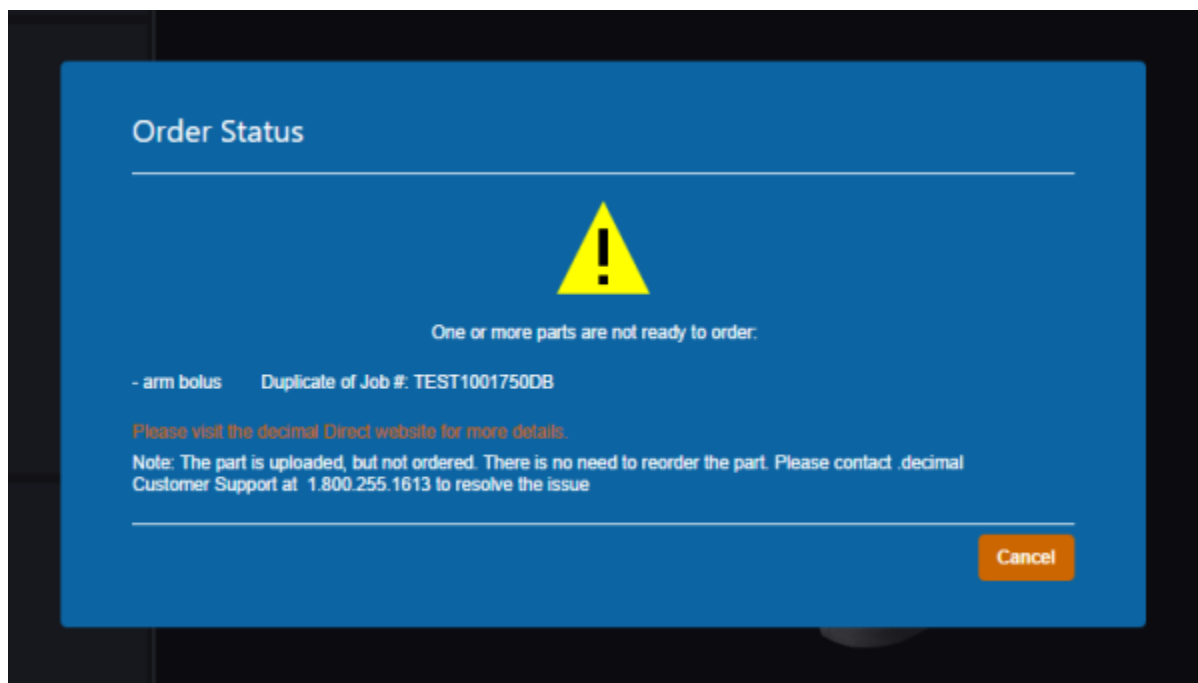


Users should contact .decimal Customer Support (1-800-255-1613) for assistance in placing their order to ensure a timely delivery.

Note: Depending on the error your order may successfully be uploaded to decimal Direct. So there may be no need to re-upload the order.

## Application Settings

The following settings are present within the decimal Ordering Tool settings UI accessed within the View → Settings menu:

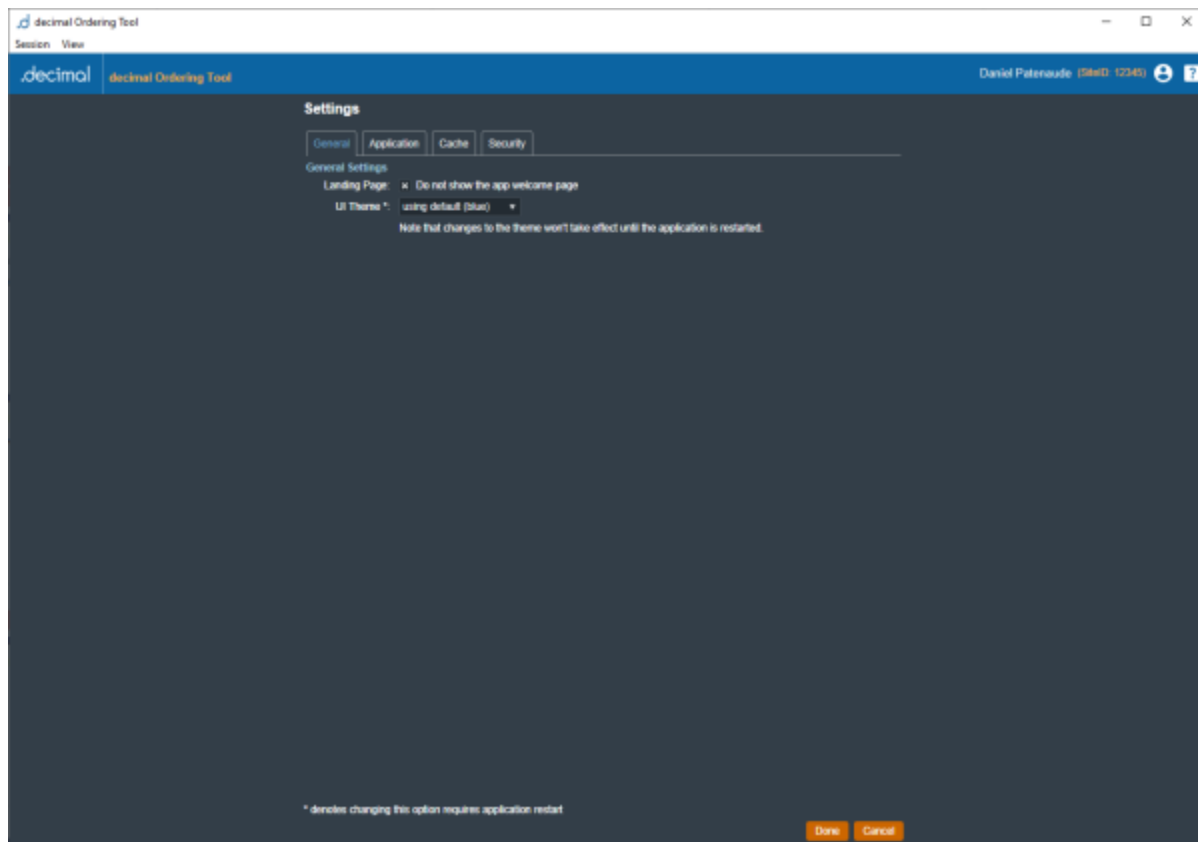

### General

- **Landing Page**
  - **Do not show the app welcome page:** Allows the user to specify whether or not to view the application landing page when decimal eRT is opened.
- **UI Theme**
  - Set the visual theme of the application (Note: the 'blue' theme is the currently supported theme; a 'dark' theme is available, but is a deprecated theme and not routinely updated at this time).

## Application

There are no specific application settings for the decimal Ordering Tool at this time.

## Cache

The local disk cache is where calculation results are stored on the client computer.

Users can choose to manually clear the disk cache or let it clear oldest data as the cache fills. Refer to [Data Management Storage and Caching](#) for more information.

## Security

- **Inactivity Timeout**
  - Specifies the time out period (in minutes) to automatically log the user out of the application. Timeout is based on mouse/keyboard interaction within the user interface. Note: This inactivity timeout is different than the user session token received from the [decimal Launcher](#).
- **SSL Settings**
  - **Disable SSL Certificate Revocation Check:** Enabling this option will disable SSL Certificate Revocation checking for proxy networks. This is not recommended unless your facility network is using a proxy that is causing SSL certificate validation errors within the application.

From:  
<https://apps.dotdecimal.com/> - **decimal App Documentation**

Permanent link:  
[https://apps.dotdecimal.com/doku.php?id=decimal\\_ordering\\_tool:userguide:tutorials&rev=1654710415](https://apps.dotdecimal.com/doku.php?id=decimal_ordering_tool:userguide:tutorials&rev=1654710415)

Last update: **2022/06/08 17:46**