

# DICOM Patient Import

In order to begin planning in the decimal ElectronRT App you must first have patient imaging captured and structure geometries defined. This information is brought in the ElectronRT App by importing DICOM CT and Structure Set files.

**Note: DICOM patient import will be disabled until the organization configuration has been completed.**

## Default Import Directory

The default import directory can be set in the app settings by navigating to View → Settings. This directory will always be the default directory when importing a patient. Users can then further refine the import folder as needed starting from this default level.

## Importing a New Patient

### Local File Import

If you have a set of Ct Images and structures you can import them directly into a new patient through the “Import Patients” menu.

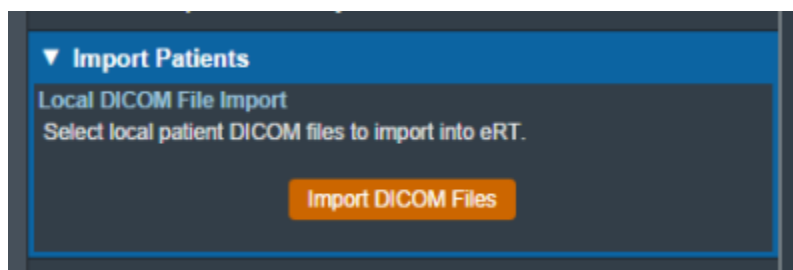


Fig. 1: Import in the Main App Page

Selecting the “Import New Patient” option will open the Importing UI.

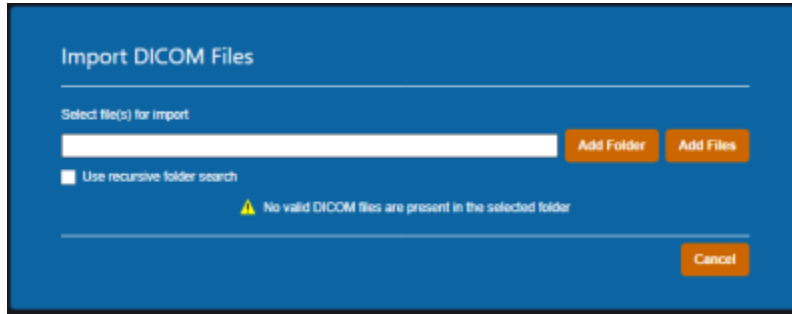


Fig. 2: Import Directory Input

Here you may either browse your computer or copy the path to the folder that contains the CT images

and structures you wish to import.



: 0.0.3 Add in new DICOM browser options.

Refer to the [Processing Imports](#) section once files have been selected to import.

### DICOM Receiver Imports



### Processing Imports

Then the import will process:

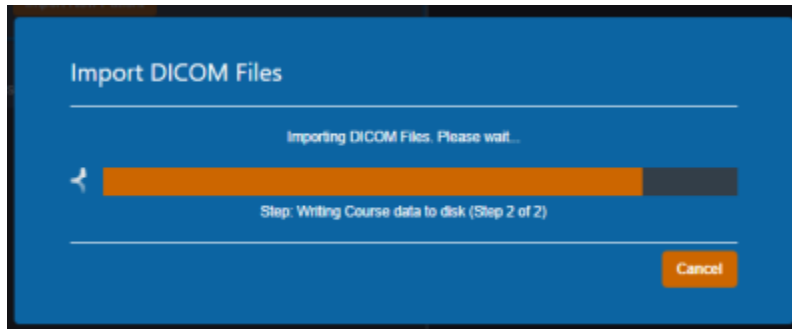


Fig. 4: Import Loading screen

If there are no Errors while importing you will be taken to a confirmation page for your import.

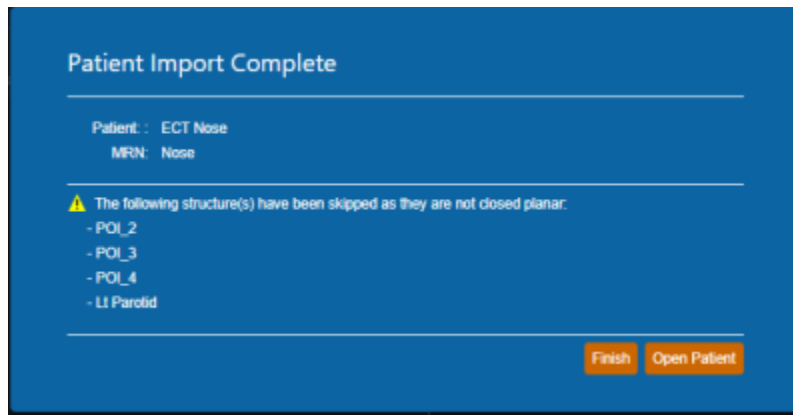


Fig. 5: Import Finalized

If there are errors please refer to the [Common Errors](#) section below.

## Common Errors

### Missing/Undefined External Structure

If the imported DICOM Structure Set does not have a structure flagged as *RT ROI Interpreted Type (3006,00A4) EXTERNAL* you will be directed to specify the external patient structure before the import resumes. The patient external structure must be set to the defining boundary of the patient's outer surface, and not a rind skin structure.

Note: You will be warned if the selected external structure is not the structure with the largest volume. This is to prevent incorrect structure selection. For example, when a 'External' and 'Skin' structure both exist, and the 'Skin' structure is a rind, if the user selects the 'Skin' structure, the warning will state the 'External' structure has a larger volume, since that's the correct representation of the patient structure needed for the ElectronRT App.

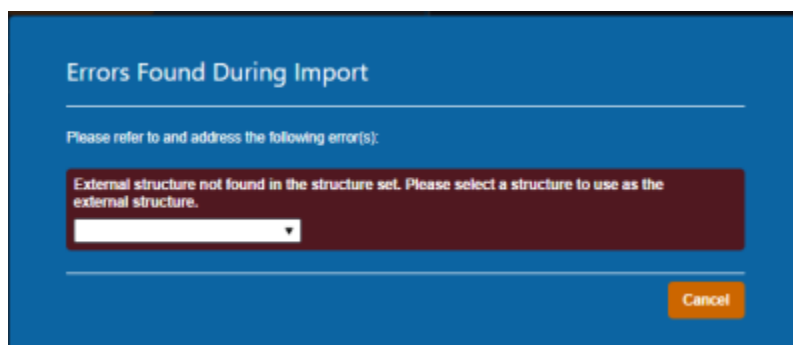


Fig. 6: Missing External Error

The drop down menu will have a list of the structures in the imported set. You will be able to select one as the external for this course then confirm your choice by pressing "Set External".

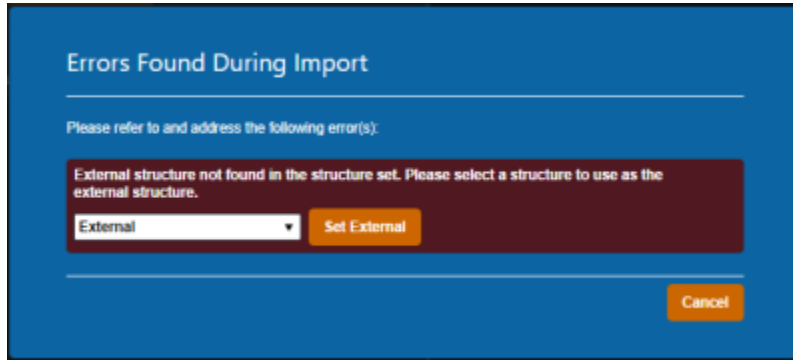


Fig. 7: Set External Dialog

## Structures Skipped During Import

Structures can be skipped during Structure Set import if the structures meet the following criteria:

- Non closed planar (e.g.: points or structures where slices are not fully closed)

Structures that are skipped will be denoted at the end of the DICOM import as shown in [figure 8](#).

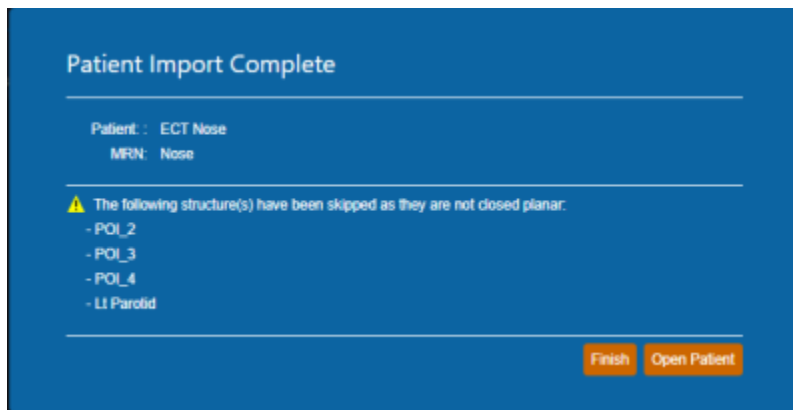


Fig. 8: Import Finalized

## Importing an Existing Patient

If the files you have selected to import coincide with the MRN of a patient that already exists in the app you will receive this error.

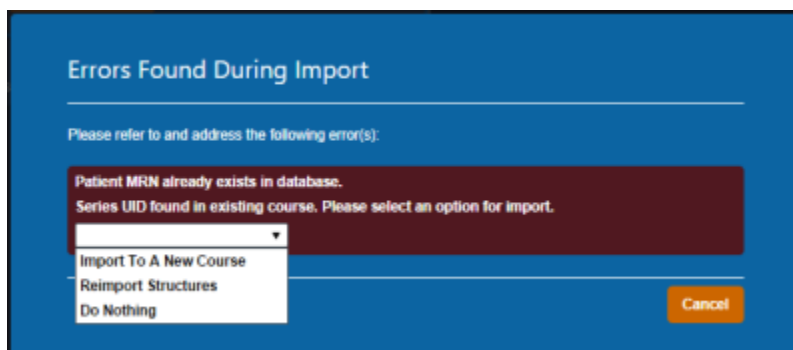


Fig. 9: Existing Patient Error

You have three options provided to resolve the error:

- **Import to a New Course**
- **Re-Import Structures**
- **Do Nothing**

### Import to a New Course

Importing to a new course will simply complete the import with the app creating a new course for this patient using the new imported files (note this option will not affect any existing courses or plans for the patient).

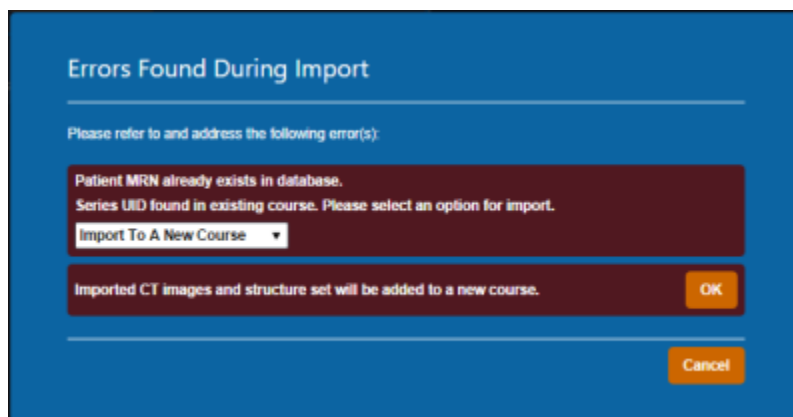


Fig. 10: Importing as a New Course

### Re-import Structures

Re-Importing structures can be used if your goal is to update or add structures to an existing patient Course.

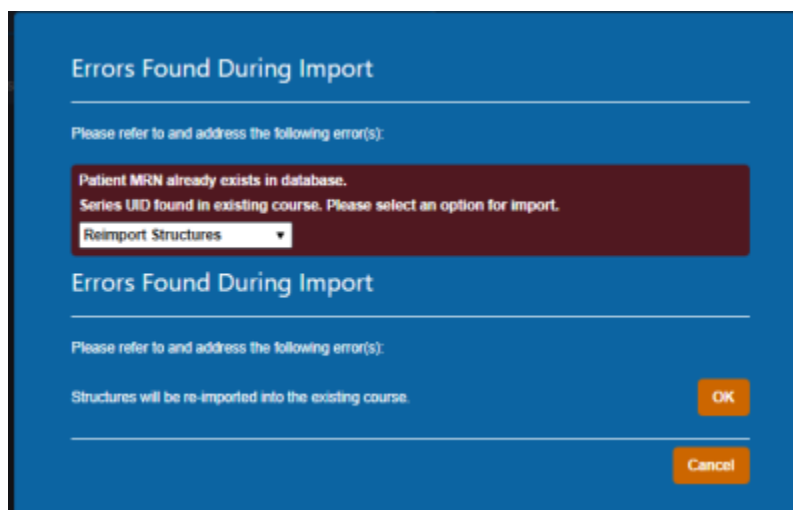


Fig. 11: Re-Import Structures

After selecting this option, you will be presented with a list of the structures in the import.

For each you may select one of three options:


- **Do not Re-import:** Skips this structure for the re import.
- **Replace existing geometry:** Replaces the existing structure in the current Course with the imported one.
- **Import as new structure:** Imports the structure to the Course as a new structure while not affecting the existing one and allows the user to specify a new name for the new structure. If the structure does not exist currently it creates it as normal.

The screenshot shows a blue dialog box titled "Errors Found During Import". It contains a message: "Please refer to and address the following error(s):" followed by a red box with the text: "Patient MRN already exists in database. Series UID found in existing course. Please select an option for import." Below this is a dropdown menu currently set to "Reimport Structures". The dialog then lists several structures with corresponding dropdown menus, all currently set to "Do Not Reimport". The structures listed are: PTV, PTV 6840, Rt Parotid, Cord, Cord+5, Ant avoid, Post avoid, Rt Eye, Lt Eye, External, and ROI US 1CM. A "Cancel" button is located at the bottom right of the dialog.

Structure Name	Selected Option
PTV	Do Not Reimport
PTV 6840	Do Not Reimport
Rt Parotid	Do Not Reimport
Cord	Do Not Reimport
Cord+5	Do Not Reimport
Ant avoid	Do Not Reimport
Post avoid	Do Not Reimport
Rt Eye	Do Not Reimport
Lt Eye	Do Not Reimport
External	Do Not Reimport
ROI US 1CM	Do Not Reimport

Fig. 12: Re-Importing Structure List

Once you have made your decision for each structure you must select "Re-Import Structures" to finalize your changes.

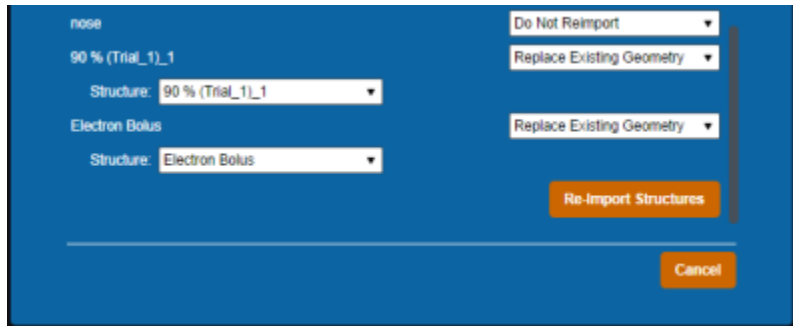


Structure	Action
nose	Do Not Reimport
90 % (Trial_1)_1	Replace Existing Geometry
Electron Bolus	Replace Existing Geometry

Buttons: Re-Import Structures, Cancel

Fig. 13: Finalize Re-Import

## Do Nothing

This is the simplest action, selecting “Do Nothing” will cancel your import and return you to the main menu.

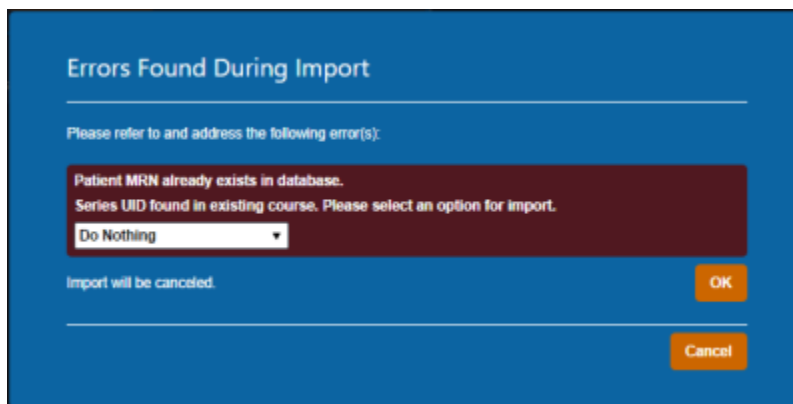

Errors Found During Import

Please refer to and address the following error(s):

- Patient MRN already exists in database.
- Series UID found in existing course. Please select an option for import.

Do Nothing

Import will be canceled.

Buttons: OK, Cancel

Fig. 14: Do Nothing

From:  
<http://apps.dotdecimal.com/> - **decimal App Documentation**

Permanent link:  
[http://apps.dotdecimal.com/doku.php?id=electronrt:userguide:tutorials:dicom\\_import&rev=1615997615](http://apps.dotdecimal.com/doku.php?id=electronrt:userguide:tutorials:dicom_import&rev=1615997615)

Last update: **2021/07/29 18:24**

