




This page is deprecated and now included in the [Astroid Optimization \(PBS\)](#) page

Optimization Objectives

Objectives communicate to the optimizer the goals that are important to strive for in your plan. *Objectives* are set at the *Plan* level under *Plan Constraints/Objectives* and they apply to the total, combined dose from all beams. *Objectives* are not given any relative importance at this point (i.e. their order within the list is not meaningful). The *Objectives* drive the solution of the Multi Criteria Optimization (MCO) and for each *Objective*, a corresponding *Navigation Slider* will be presented to allow for exploration of trade-offs in the case of competing objectives (for more information about the MCO process and how objective importance/weighting is handled in Astroid refer to [this article](#)).

The following objective selections are available in Astroid:

- **min_max**: Minimize the maximum dose within a structure (drive dose down)
- **max_min**: Maximize the minimum dose within a structure (drive dose up)

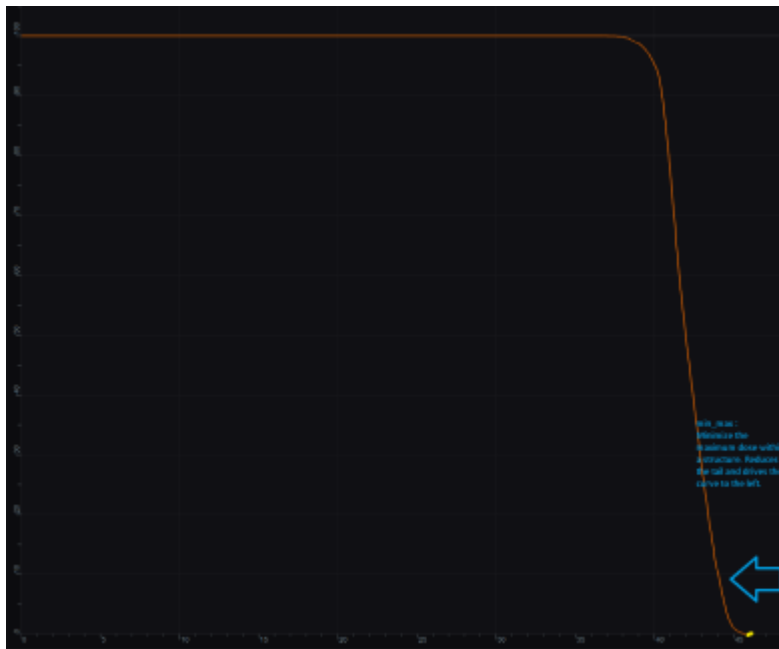

min_max: Minimize the Max Dose




max_min: Maximize the Min Dose

- **min_mean**: Minimize the mean dose within a structure (drive dose down)
- **max_mean**: Maximize the mean dose across the structure (drive dose up)
- **min_overdose**: Minimize the high dose within a structure
 - Dose will be driven down only until the specified limit is reached (this is often more relevant that min_max, since it may not be beneficial to continue minimizing beyond a certain dose level)
- **min_underdose**: Minimize the low dose within a structure
 - Dose will be driven up only until the specified limit is reached (this is often more relevant that max_min, since it may not be beneficial to continue maximizing beyond a certain dose level)



min_overdose: Minimize the high dose

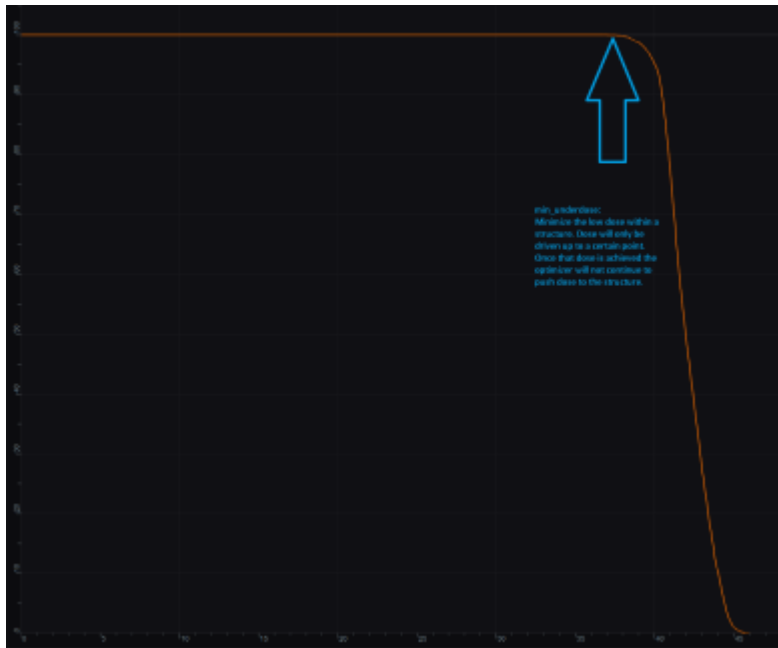

 min_underdose: Minimize the low dose

Working with Objectives

1. Open the *Objectives/Optimizer* sub-block contained in the *Optimization* block
2. Choose a structure to which you wish to apply objectives
3. Check the boxes to activate the desired objectives for the structure and then set the dose level if applicable

Objectives

☒ Rectum ☒

☒ min_max

☐ min_mean

☒ min_overdose

36 Gy(RBE)

☒ PTV ☒

☐ min_max

☐ min_mean

☒ min_overdose

46 Gy(RBE)

☐ max_min

☐ max_mean

☒ min_underdose

52 Gy(RBE)

Add Structure ▼

OK Cancel

Once all the *Objectives* have been set, the user is ready to run the MCO solver, which is performed in the *Objectives/Optimizer* block.

From:

<https://apps.dotdecimal.com/> - **decimal App Documentation**

Permanent link:

https://apps.dotdecimal.com/doku.php?id=planning:userguide:tutorials:optimization_objectives

Last update: **2021/07/29 18:28**