


Creating a PBS Beam

Defining treatment beams will be one of the most important tasks within the Astroid planning system. Defining appropriate beams will require users to use their knowledge and experience to properly select many of the parameters that define a treatment beam. These parameters include the target, geometry (isocenter, gantry and couch angles), beamline devices, air gap, and spot placement options. The *Beam* task utilizes a series of blocks to organize the beam creation process into a common step-by-step sequence. Several blocks are optional as not all beams will use all features. Additionally, it is important to point out that the treatment room & default spot placement parameters are set outside of the individual beam creation tasks as these apply to all beams (however, spot placement parameters can be overridden within each beam if desired). An example of constructing a lateral beam, with the isocenter at the centroid of the PTV is given below to illustrate the features available when defining a beam.


1. From within the *Plan Overview* select the *Beams* block




2. Select the treatment room from the drop down then click the blue OK button

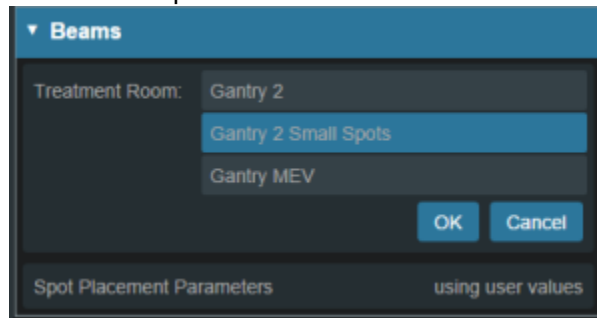

3. Next select the *Spot Placement Parameters*. Edit these parameters to define the spot placement grid for each beam in the plan, noting that individual beams can override these values during beam creation if you so desire

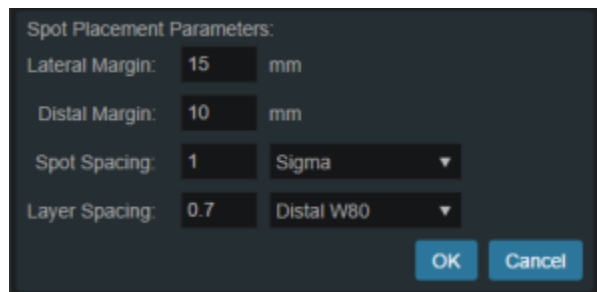

4. The next step is to create the beam by clicking the *Create New PBS Beam* button

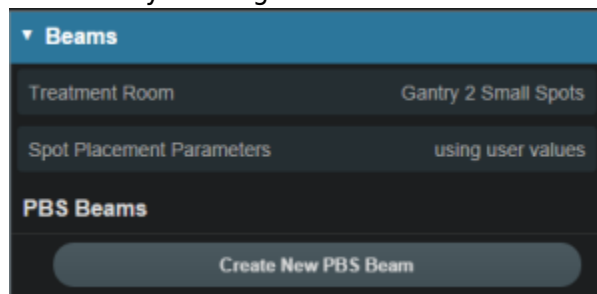


5. Now we will proceed step-by-step through the various “blocks” to create a complete beam as shown below:

1. In the *General* block select the intended target from the dropdown. You may choose an existing target or create a new structure. For this example we chose the PTV7920 as the target. The target is important because the target will be used to define the extents of the spot placements as well as the aperture (if used)




2. In the *Approach* block specify the desired isocenter from the dropdown. You may choose to use the centroid of the target (as shown below) or you can select or create a new point to define the location for the isocenter. You will also enter the gantry angle and couch angle here as well. Editing these values can be done by typing directly in the provided fields or by using the sliders. The patient in this example is feet first so we will use 90 for the *Gantry* angle and 0 for the *Couch* Angle to create a left lateral beam

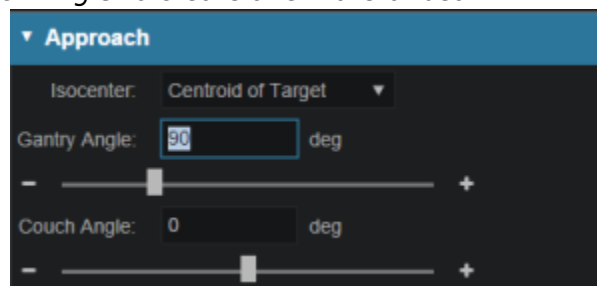

3. The next block is the *Snout* block. A list of snouts associated with your specified treatment room will be available to choose from. The 12cm snout was chosen here.



4. If desired, add an aperture in the *Aperture* block. (No *Aperture* was chosen in this example)
 1. Refer to [Creating a New Aperture](#) for detailed instruction
5. If desired, select the range *Shifter* to use based on the ones available for the selected snout. If no *Shifter* is needed as in the example given, the user may go to the next step

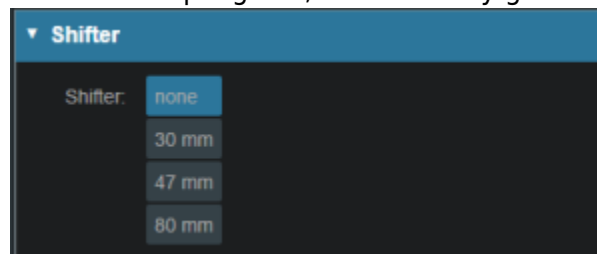

6. Once the beam line devices are defined, we can move to specify the *Air Gap* distance. The valid air gap range will be listed based on the selected snout. The user may choose any value in this range. 30mm was chosen in the below example.

Create PBS Beam

▼ General

Color: [11 colored circles]

Label: G0 C0; no snout; no shifter;

☒ automatically generate label

Target: [dropdown]

Description: [text area]

► Approach none

► Snout none

► Aperture none

► Shifter none

► Air Gap 30 mm

► Spot Placement using plan settings

► Proton DRR Options

7. With the beam positioned and any beamline devices put in place, the user is ready view the PBS Spots and adjust the *Spot Placement* values if needed. The *Spot Placement* box, if chosen, will allow the user to set new parameters, overriding the spot placement parameters for this one beam if desired. The example below illustrates the message shown when using the spot placement values from the plan level.

▼ Spot Placement

Plan spot placement settings being used.

[Use beam spot placement settings](#)

8. The last block in the *Create Beam* task is the *Proton DRR Options*. These do not impact the beam and are used purely for visualization purposes so that you can set the DRR Options to levels that generate appropriate anatomy visualizations. An example DRR is shown below. Note that Astroid allows you to define 2 distinct DRRs and then blend them together using a simple weight factor to create a single DRR image on the screen. This gives users the freedom to create high contrast, high quality DRR visualizations.

▼ Proton DRR Options

Window: 1500 Level: 750

Hu Settings 1

Weight: 1

Bone ▼

Min: 200 Max: 2995

Advanced Options

Hu Settings 2

Weight: 0

Soft Tissue ▼

Min: 40 Max: 80

Advanced Options

Image Z Position: -200 mm

From:
<https://apps.dotdecimal.com/> - **decimal App Documentation**

Permanent link:
https://apps.dotdecimal.com/doku.php?id=planning:userguide:tutorials:pbs_beams&rev=1487010231

Last update: **2021/07/29 18:25**