

Beam Delivery (SOBP)


SOBP Fraction Groups

Defining *Fraction Groups* is the final step in the SOBP planning process within Astroid. A fraction group is simply an arrangement of beams that will be used in a typical daily treatment fraction. The Fraction Group contains a list of beams where each beam can be weighted and normalized based on the needs of the patient at hand. Beam normalizations allow for minor field independent adjustments of dose and these factors are limited to stay within 0.90 - 1.10. Beam Weights established the relative fluence of each beam within the fraction group and therefore the sum of the beam weights in each fraction group must be 1.0. Further details of defining fraction groups are provided below.


- To get started, navigate to the *Beam Delivery* block and click *Create New Fraction Group*.




General Fraction Group Data

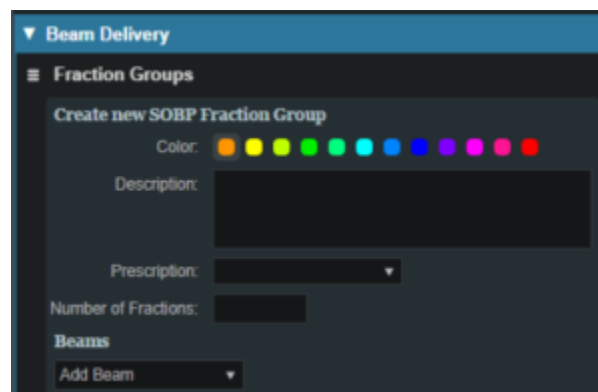

- **Color:** Display color of the *Fraction Group*
- **Description:** Optional, user specified text describing the *Fraction Group*
- **Prescription:** Prescription that the *Fraction Group* implements; note that only targets containing dose statements from this Prescription will be available when selecting the Target for the Fraction Group

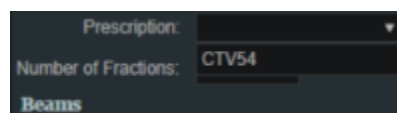


- **Number of Fractions:** The total number of fractions to be delivered for this *Fraction Group*; this is very important as it will determine the appropriate Monitor Units for the individual beams

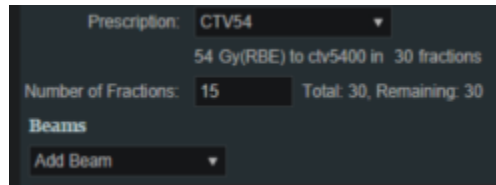

Beam Options

- Select the beams to be added to this fraction group from the drop down list. Repeat this process as necessary. Notice how a one the beams in the image below is grayed out. This means that the beam has already been selected or it is a *patch beam*. If it is a patch beam, the beam will automatically be added with its parent beam. There is also the option to create a new beam from the drop-down list.

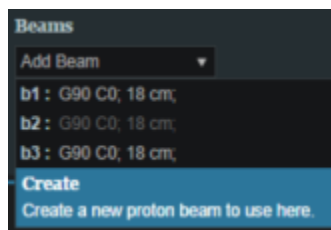

- Change the normalization factor by clicking its box and entering the value desired.
 - This is independent from the weight values.
 - Available as an option for patch beams
 - Useful when adjusting hot or cold areas where beams treat the same area or meet.
- Weight each beam by either sliding the slider or clicking on the weight value box and type in the value.
 - The value of all beams listed **MUST add up to 1.0** (100%).
 - By default, the beams are automatically weighted evenly, and, if edited, can be reset by clicking the *weight beams evenly* option.
 - Use the *Beam Lock* (not available for patch beams) to restrict a beam from automatically adjusting when another beam's weight is adjusted. This option is useful to secure one beam's weight while adjusting the other beams' weights (see image below). The sum of all weights still must equal zero.

Beam Delivery

Fraction Groups

Create new SOBP Fraction Group

Color: [Color Selection]

Description: [Text Field]

Prescription: CTV54
54 Gy(RBE) to chv5400 in 30 fractions

Number of Fractions: 15 Total: 30, Remaining: 30

Beams

b1 : G90 C0; 18 cm; Normalization Factor: 1

Beam Weight: 0.5

b2 : G90 C0; 18 cm; (patch) Normalization Factor: 1

b3 : G90 C0; 18 cm; Normalization Factor: 1

Beam Weight: 0.2

b4 : G90 C0; 18 cm; Normalization Factor: 1

Beam Weight: 0.3

Add Beam [Dropdown]

weight beams evenly

Done Cancel

- When finished, click the *Done* button. Clicking the *Cancel* button will reset the changes made to the current fraction group and return the *Beam Delivery* block to its previous state.
- Repeat process of adding fraction groups until all remaining fractions are completed.

From:
<https://apps.dotdecimal.com/> - **decimal App Documentation**

Permanent link:
https://apps.dotdecimal.com/doku.php?id=planning:userguide:tutorials:sobp_beam_delivery&rev=1569005407

Last update: **2021/07/29 18:25**