2025/12/15 17:42 1/5 IT/Network Troubleshooting

IT/Network Troubleshooting



.decimal System Status

Prior to beginning the network troubleshooting steps refer to our status page to determine if there is a service outage with any of the .decimal backend services.

As the software provided by .decimal requires access to the internet to place orders, the proper functioning of our software can be affected by the local network within which the software is installed. Since each customer facility implements their own unique security, firewall, and networking policies we cannot guarantee support for all configurations. Therefore, this guide was created to provide general troubleshooting steps to help our customers and their IT departments ensure they are not blocking or interfering with .decimal software and the HTTPS requests to .decimal's servers that are required to install, authenticate, and place device orders.


Traffic Architecture

.decimal's applications communicate via HTTPS to our server resources. The diagram below demonstrates how our HTTPS traffic flows through each piece of infrastructure and identifies common failure points:

2025/12/15 17:42 2/5 IT/Network Troubleshooting



.decimal Traffic Architecture

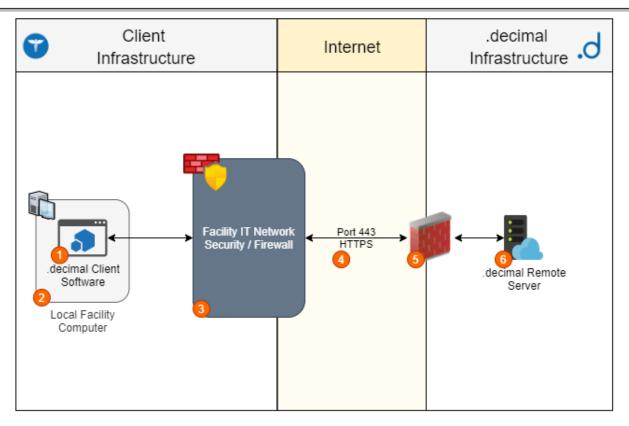

Notes

- decimal client software is provided by .decimal and installed on your network computers.
- 2 The customer's local PC in which the .decimal client software is installed. This may include local Antivirus or software firewalls that can hinder and impede .decimal software and traffic if not configured correctly.
- 3 The customer's unique security and firewall configuration. This may include any number of network firewalls, web security appliances, proxies, etc. .decimal client software traffic passes through your facilities IT security devices, which can impede and hinder that HTTPS communication. Each device within your facility's network security suite that .decimal's traffic flows through should have an exception for each required .decimal url/port.
- 4 .decimal HTTPS traffic can have multiple destinations depending on the software being installed. Refer to the software's Network and System Requirements document provided by .decimal Staff. At a minimum the URL direct.dotdecimal.com (64.128.252.104) port 443 is always used.
- .decimal's internet firewall. We do not blacklist or filter our own HTTPS traffic.
- 6 .decimal's server infrastructure within our network. Our servers accept your HTTPS traffic as received from .decimal's client software after it exits your network security devices.

Common Issues and Resolutions

1. Blocking or filtering the required .decimal urls or ports:

.decimal applications require specific urls(or IPs) and ports to be unblocked and unimpeded

2025/12/15 17:42 3/5 IT/Network Troubleshooting

in order to authenticate and place orders to .decimal's servers. Refer to the System and Network Requirements document provided by .decimal staff. At a minimum the url direct.dotdecimal.com (64.128.252.104) port 443 is always used.

2. Proxies / Web Security Appliance (WSA):

 Some WSAs will attempt to decrypt and inspect our HTTPS traffic. Some customers have needed to include an exception in this network security device as this can impact and cause our encrypted and authenticated HTTPS traffic to encounter errors while in transit between the client software installed at the clinical facility and .decimal's servers.

Troubleshooting Steps

The following steps may help when troubleshooting network security and IT issues for your facility.

1. Compatibility Checker


Running the compatibility checker will ensure that all required .decimal IP addresses and ports are not being blocked.

Note:



The Compatibility Checker only ensures that the primary .decimal IPs and ports are not being completely blocked or blacklisted. Additional network security devices may impede or be filtering the HTTPS traffic. If .decimal software continues to encounter problems, please ensure other network security appliances and devices are not filtering or impeding .decimal's HTTPS traffic.

2025/12/15 17:42 4/5 IT/Network Troubleshooting

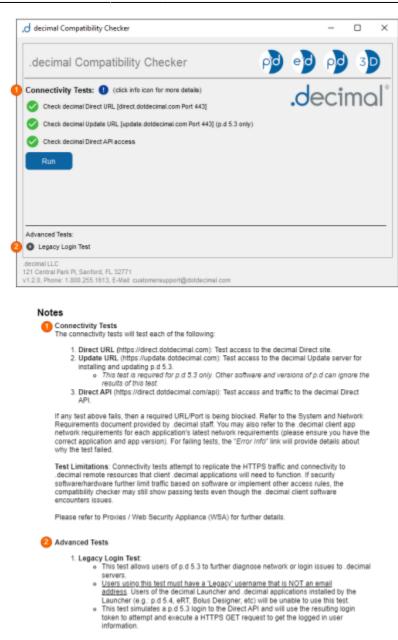

2. Guest/Unrestricted Network

This option is advised if your network IT has whitelisted and unblocked .decimal's urls and ports but you still encounter issues using .decimal's software and communicating to .decimal's servers.

- 1. Attempt to download and install the .decimal client software on a device (e.g.: a laptop) not connected to your main cooperate network (e.g.: a guest WiFi network or smart phone hotspot) that has absolutely zero security or firewall blocking, filtering, or packet inspection.
- 2. Attempt to login to the .decimal client software
- 3. Confirm you are able to download, install, login, and use the .decimal software on an unrestricted device on an unrestricted internet connection. This ensures no security, firewall, or packet inspection security policies are interfering with your HTTPS connection and traffic to .decimal servers.
 - 1. **Conclusion:** If Step #3 passes, then there is still security policies in place on your main

2025/12/15 17:42 5/5 IT/Network Troubleshooting

facility network impeding .decimal HTTPS traffic.

From:

http://apps.dotdecimal.com/ - decimal App Documentation

Permanent link:

http://apps.dotdecimal.com/doku.php?id=support:it_troubleshooting&rev=1637241521

Last update: 2021/11/18 13:18

